



User's Manual

24-Minute Miracle • Bright Skin for 24 hours!



TABLE OF CONTENTS

Product Introduction.....	4
Product Principle	5
Product Composition.....	7
User’s Manual [details].....	8
User’s Manual [operation manual].....	10
Applications (Dermatology Clinic)	12
Certification	13
Warranty.....	14

With the 4 sources of the light effect of OzMask
and synergistic effect of **light and heat!**

IMPROVE SKIN WITH LIGHT!



*Improving freckles, skin problem, skin elasticity,
wrinkle, and aging process set free women problems!*

Thank you for purchasing beauty rock.


To ensure a more efficient and safe use of this product, please read and understand the precautions before using it.

PRODUCT INTRODUCTION

This product is a new concept of PDTT(Photo Dynamic Thermo Therapy) esthetic device that combines PDT(Photo Dynamic Therapy) and Thermo Therapy

It is a new-concept skin beauty device which is a collaboration of the detox and the effects of cell regeneration of far infrared rays in addition to the healing powers for each color wavelength of BRY-G using LEDs.

PDT
[Photo Dynamic Therapy]
+
Thermo Therapy
||
PDTT
[Photo Dynamic Thermo
Therapy]

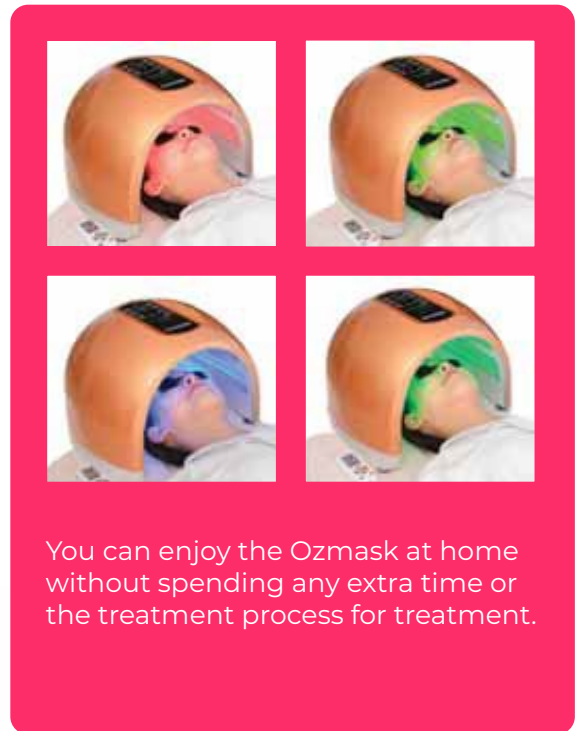

PRODUCT PRINCIPLE

24-Minute Miracle • Bright Skin for 24 hours!

LED(RGBY) for LLLT(Low Level Laser Therapy) is a very stable light source that treats skin diseases and problems by varying indications depending on the wavelength(nm) of the light source without damaging the skin with a warm light.



BLUE 455nm | Depth Of Penetration ~ 1 mm

- Inhibition of sebum production • Prevention of wound infection
- Sterilization of bacteria in acne • Calming skin problems



GREEN 530nm | Depth Of Penetration ~ 1 mm

Helps in elastin production

- Improving skin problems • Calming skin problems



Yellow 580nm | Depth Of Penetration ~6mm

Helps in collagen and cell reactivation

- Skin capillary contracture • Improving facial blushing • Collagen production in skin • Improving skin elasticity • Improving skin wrinkles



RED 630nm | Depth Of Penetration ~12mm

By activating enzymes in underlying skin layers, It activates mitochondrial metabolism and helps produce collagen and elastin, helping with cell regeneration and reactivation

- Skin cell activation and reactivation • Improving skin elasticity • Improving skin wrinkles • Skin cell regeneration • Blood circulation • Relieve acne scars




FIR-RAY 970mm~4µm | Depth Of Penetration ~40mm

By transferring heat energy to the subdermal blood vessels and fat cell layers, It improves blood flow and induces transparent skin by helping to eliminate wastes and toxins in fat cells.

- Improving blood flow • Improving skin tone • Improving skin elasticity • Skin detoxing

The effects of OzMask which improve and prevent skin aging Penetrating Mimics of Healing Rays

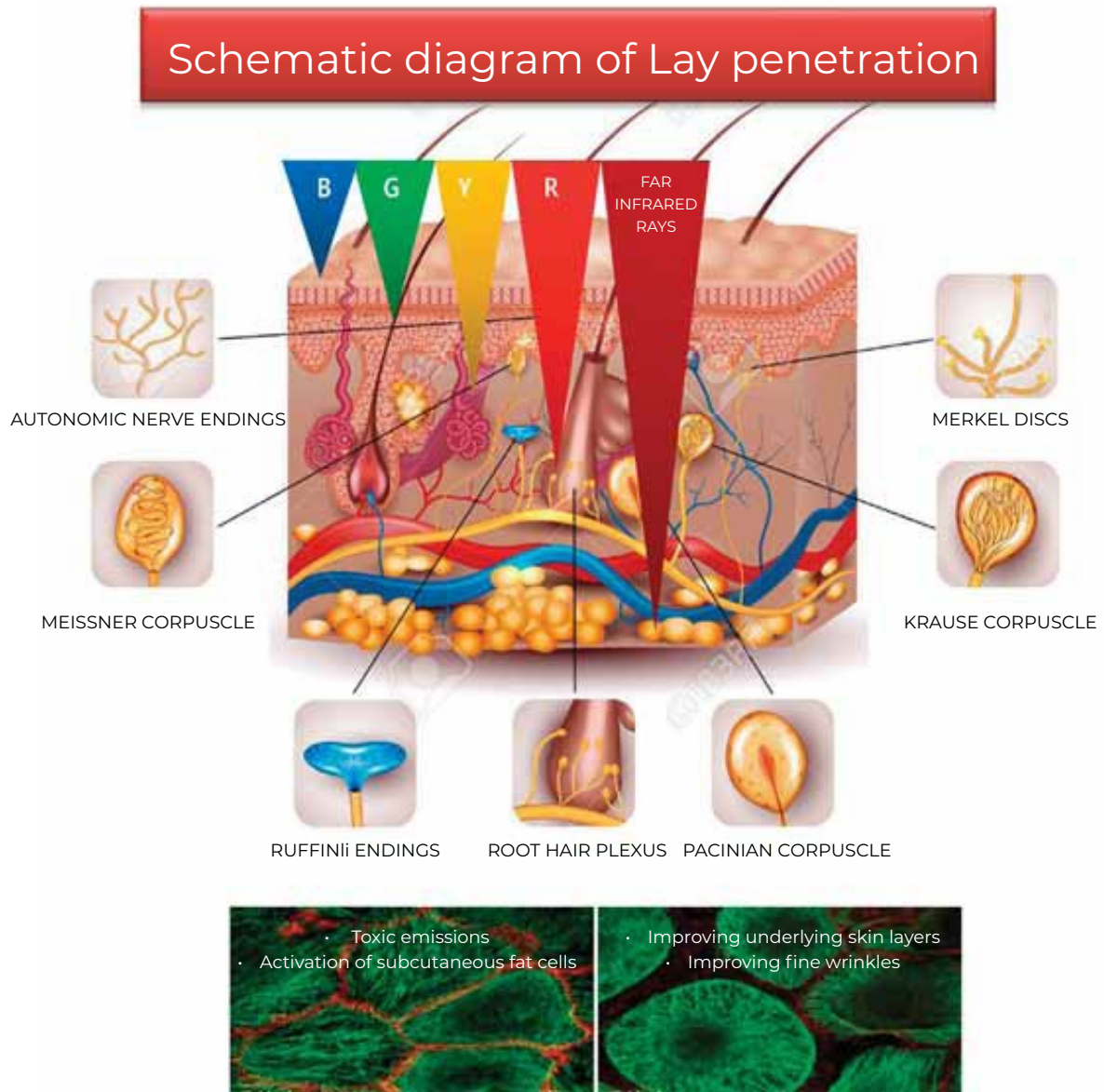

BLUE

effects of acne treatment, relieving atopic dermatitis

RED

relieving wrinkles, skin moisturization, strengthening skin elasticity

YELLOW

collagen production, discolored skin treatment with light, Improving skin resistance against UV

GREEN

relieving atopic dermatitis, suppression of fungi, prevention of food poisoning

FIR-RAY

cell growth promotion, wound healing, pain treatment, muscle relaxation

PRODUCT COMPOSITION

OzMask Main Body



Adapter with Power Cable



DC 24V / 5A

Tencel Mask pack & Mist



(1 BOX/15pcs)

This product can be repurchased separately after use

User's Manual




Eye Mask (Silicon)



(2pcs)

PRODUCT STRUCTURE

1. Non-slip rubber packing
2. DC Power Inlet(Adapter)
3. Far Infrared Heater
4. RYGB Focusing Lens product

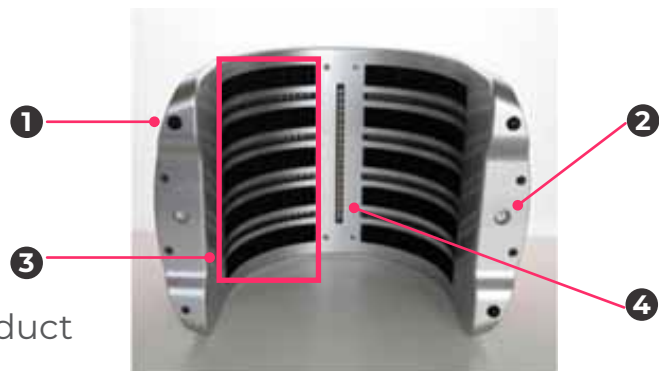

bottom Sight


USER'S MANUAL (DETAILS)



Active part viewed from above the product

If you remove the protective film vinyl attached to the plaque of the active part at the time of use after the first purchase, the product will be much more luxurious

Auto Mode

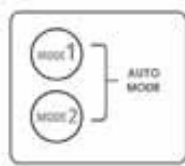

If you press the text button written as **Mode 1**, it works without pressing the ON button.

Mode 2 is the same. Only slightly touch the ON/OFF button to turn it off.



Auto Mode 1

Slightly press button **Mode 1** of auto mode at the time of use.

Auto Mode 1 operates in the order of red, yellow, green, and blue (RYGB) color for **about 24 minutes**. The initial temperature changes from 70°C to 60° C 10 minutes after an operation and there will be sound by changes.

Do not be surprised because this is not a warning horn for any trouble or operational failure.

Auto Mode 2

Slightly press button **Mode 2** of auto mode at the time of use.

Auto Mode 2 operates in the order of red, yellow, green, and blue (RYGB) color for **about 35 minutes**. The initial temperature changes from 70°C to 60°C 10 minutes after an operation and there will be sound by changes.

Do not be surprised because this is not a warning horn for any trouble or operational failure.

The difference between **Auto Mode 1** and **Auto Mode 2** is only 24 minutes for **Auto Mode 1** and 35 minutes for automatic mode 2.

There is only a difference in operating time. Consumers can use it conveniently by themselves.

We recommend you use auto mode only!

Recommended Mode: Auto Mode

Heater Manual Mode

- LED and thermal mode can be used at the same time and the illumination intensity can be adjusted in order of low/medium/high level by touching the (RYGB).
- When touching the ON button, the initial set-up temperature is 70° C and the operating time is set to 60 minutes.
- The temperature is adjusted in 1 degree and it is adjusted by 1 degree when pressing continuously.
- The time is adjusted in 5 minutes and changes by 5 minutes when pressing continuously

LED Manual Mode



If you press the button as shown in the picture, you can adjust the temperature step by step, but we recommend you only use auto mode before you are used to using the device if possible.

Heater Manual Mode



- LED and thermal mode can be used at the same time
- The illumination intensity can be adjusted in order of low/medium/high '0: level by touching the (RYGB)

USER'S MANUAL (OPERATION MANUAL)

Thank you for purchasing beauty rock.

To ensure a more efficient and safe use of this product, please read and understand the precautions before using it.

1. Power Operation: Connect the product and the power cable.

- Power-Supply Specification: DC 24V/5A.
(Use the specified power supply. If not, it could cause a failure)
- "--" mark appears in the temperature setting window and the time setting window when the initial power-on.

2. The setting of LED operating mode

1) AUTO MODE

- MODE 1 operates in the order of RYGB(red, yellow, green, and blue) for 24 minutes. The initial heating temperature is 70°C and the time will be displayed.
- MODE 2 operates in the order of RYGB(red, yellow, green, and blue) for 35 minutes. The initial heating temperature is 70°C and the time will be displayed.
- Touch auto mode or ON/OFF to turn off the auto mode.

2) MANUAL MODE

- LED and thermal mode can be used at the same time
- The illumination intensity can be adjusted in order of low/medium/high level by touching each button of the (RYGB).

3. The setting of temperature/time in controlling heating mode

- When touching the ON button, the initial set-up temperature is 70°C and the operating time is set to 60 minutes.
- The temperature is adjusted in 1 degree and it is adjusted by 1 degree when pressing continuously.
- The time is adjusted in 5 minutes and changes by 5 minutes when pressing continuously.

4. AUTO MODE 1

- Pressing the Mode Key initiates operation and displays the current temperature in the temperature display window.
- The far infrared heater operates at a maximum of 70°C and at the same time it flashes every 10 seconds of RED, 8 seconds of YELLOW, 7 seconds of GREEN, and 5 seconds of BLUE.
- The set-up temperature is automatically changed to 60°C 10 minutes after an operation and there will be beeping sound. (this is not a warning horn for any trouble or operational failure, do not be surprised)
- "MODE 1" has a total operating time of 24 minutes and is automatically turned off after the operating program.

5. AUTO MODE 2

- Pressing the Mode Key initiates operation and displays the current temperature in the temperature display window.
- The far infrared heater operates at a maximum of 70°C and at the same time it flashes every 3 seconds of RED, 3 seconds of YELLOW, 3 seconds of GREEN, and 3 seconds of BLUE.
- The set-up temperature is automatically changed to 60°C 15 minutes after an operation and there will be beeping sound. (this is not a warning horn for any trouble or operational failure, do not be surprised)
- "MODE 2" has a total operating time of 35 minutes and is automatically turned off after the operating program.

Summary of User's Manual for Operation Manual

1. Connect the product and power cable.
2. Plug the connected product's electrical cord into the power supply.
3. Slightly press button mode 1 or mode 2 of auto mode at the time of use. Heater temperature and time are automatically set.
4. Press the ON/OFF button only to turn off the power. (we recommend you only use auto mode before you are used to using the device if possible)
5. After using the product, remove the power cord from the outlet.

APPLICATIONS (DERMATOLOGY CLINIC)

Clinical Test

Evaluate abnormal symptoms that occur after using the test product



TEST
연필리버블 미백효과를 눈으로 확인하세요!

Tested by P&K Skin Research Center

Clinical test results

Increased dermal density 26% | Increased skin elasticity 13.56% | Reduced fine wrinkles 10.5% | Reduced eye wrinkles 10.36% | Reduced facial area 4.12% | Increased skin moisturization 14.85% | Improved skin tone 14.2%.

Test result

In this trial, women aged 29 to 59 years old were asked to use the test product for 4 weeks and the eye wrinkles, skin moisturization, elasticity, facial lifting, dermal density, and skin tone were measured before using the product, 2 weeks after using the product, and 4 weeks after using the product.

After the test, we obtained the following conclusions by completing the survey assessment of the test subjects.


Result

4 weeks of use of the product helps to improve skin aging (eye wrinkles, facial lifting, dermal density), skin texture, and skin tone.



CERTIFICATION

특허 출원 제: 10-2016-0112552
 디자인출원 제: 30-2016-0042993

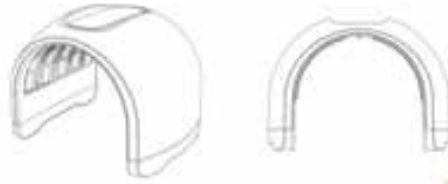

상표등록


전자파인증

P/L보험

FCC 인증

CE 인증

원격외선 시험성적서

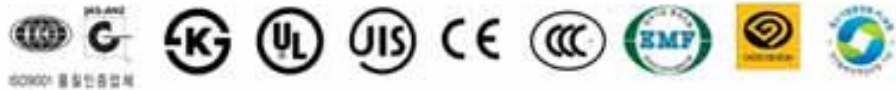

PDT TREATMENT FDA APPROVAL

PDT(Photoflin)

PDT(Levlan)

WARRANTY

Product Name	LED Skincare Device(OzMask)
Model Number	BTR 1004
Rating	DC 24V/5A
Product's Purchase Date	20 / /
Purchase Location / Phone Number	
Free warranty period	1 year from the date of purchase of the product

Please make sure to present this warranty when you receive the service.

1. The compensation standard for this product is in accordance with the compensation rule for the damage notified by the Consumer Protection Board.
2. This product is manufactured through strict quality control and inspection process.
3. We will provide free repairs if there is a failure of the product at the time of normal use during the warranty period.

In the following cases, it does not correspond to the quality condition, so it will be charged for repairs.

1. In the case of a failure due to the consumer fault or mishandling
2. In the case of a failure caused by the natural disaster (fire, earthquake, wind and flood damage, thunder stroke, etc.)
3. In the case of a product failure that occurred after the warranty period expired
4. Natural discoloration caused by the thermal function due to the nature of the product






No free repair or exchange for failure caused by reforming or modulation of the product at consumer's discretion



OzMask

24-Minute Miracle • Bright Skin for 24 hours!



 www.ozmask.com  +61 1300 000 998  OZMASK  Ozmaskau
 Penny World PTY LTD - Suite 903 10 Help Street Chatswood NSW 2067