



Professional
MACHINE

I Description

This body light-emitting diode light-weight medical aid skin beauty instrumentation has multiple gentle lights to treat totally different skin problems with different body elements. And this wonderful feature makes it a marvel skin electrical device machine. created with invasive targeted medical aid technology, this body light-emitting diode light-weight medical aid skin beauty care is best to enhance rough pores, pigment spots, to color skin, fight wrinkles, aging, skin dry, oily skin, skin allergies, skin redness, acne pits, birthmarks and fine lines. To possess the simplest results while not damaging your skin, here is a unit few must-to-read safety precautions for all the users

- Maintain a distance between a machine and the other ignitable substance around you.
- Make bound to connect the facility plug with dry hands.
- Never try and modify the instrument on your own.
- Entangle the wires if twisted before use.
- Keep the merchandise in safe, cool and dry place. Particularly avoid keeping it beneath direct daylight or in any wet place; it should injury the sunshine refulgent impact of the machine.

OzMask™ are often used alone for primary skin rejuvenation* and together with OzMask™ for wound healing, disease of the skin and for the reduction of per orbital wrinkles. Are often employed in pic Dynamic medical aid (PDT) for the treatment of non-melanoma skin cancers and severe disease of the skin. Employed in combination with OzMask™ for gentle to moderate acne and cystic disease of the skin.

I Output intensity

- 100W
- Output wavelength
- 650nm (red light)/850nm(infrared light)/415nm (Blue light)/590nm (yellow light) normal dose
- 5min
- Treatment time (standard dose) 20mins whenever
- Dose vary (adjustable) 5-60min
- Bandwidth
- what's your mean?
- Weight
- 34KG(Carton packaging)
- Unit dimensions (HxWxD) 61*42*141cm
- Head dimensions overall (Laws) 550mm*900mm
- LED head active dimensions (Laws) 880mm*270mm

I This peel tricky resolving instrumentation has numerous cultured Christmas lights including:

- RED – Wavelength of 630 nm – It the whole thing on skin reconstruction through removing furrows, aging and gloomy spots. Also, it expands the skin form, repairs the ruined skin, shrinks the pores and develops the coating snap to an exceptional extent.
- BLUE – Wavelength 470 nm – This light-weight works for those that have disease of the skin prone skin. It destructs the harmful & infected skin tissues and microorganism like- Propionibacterium.
- GREEN – Wavelength of 520 nm – together with this, it's used for equalization the oil or moisturizing components on the skin. It repairs the injury done by mental stress, dredging body fluid and evacuation.
- PURPLE – Wavelength of 690 nm – It imposingly treats the stubborn disease of the skin spots that sometimes last forever.
- CLEAR BLUE light-weight – Wavelength of 490 nm – It improves the energy of skin cells that leaves a refulgent glow on your face.
- YELLOW – Wavelength of 590 nm – It provides the skin ample quantity of element decomposes the pigments, improves the skin roughness, quality, wrinkles, redness, and moss.
- LASER light-weight – Wavelength of 510 nm – that stimulates the uptake, therefore molders the tinges, fine lines, and dusky bad skin solely.



Benefits of LED/DPL treatments

- Improves skin tone and clears blemishes.
- Uneven skin texture lowers coarseness and condenses pore extent.
- Maximizes the miens of fine lines, crumples, crow's feet, smile defenses, nasolabial bends, and anxiety lines.
- Regenerates skin by activating formative cell cells to extend production of albuminoidal and albuminoidal.
- Helps improve the looks of sun-damaged skin and irregular pigmentation.
- Enhances skin wetness that helps skin regain its immature fullness.
- Increases circulation that helps cut back animal pigment production, the explanation for brown age spots.
- Promotes nutritionary components existing among the skin.
- Stimulates and activates metabolic operate in skin cells.
- Reduces skin degradation caused by age and environmental exposure.
- Reduces the looks of overall redness and flushing.
- LED light-weight listing

So, it's best to visualize along with your doctor before you are trying it.

Product Features


We Provide ODM or OEM Service

Activate Hair Follicles to Regenerate Hair

650nm red light saves the hair follicle



Prevent hair loss and promote hair growth from the roots of the hair

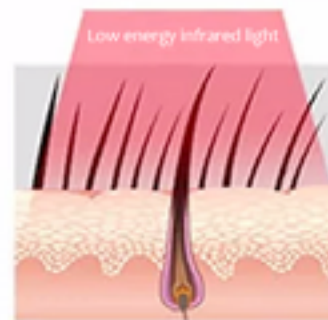

Promotes blood circulation, oil control and sewage (DHT)
improve the scalp environment for hair growth

I You are not a candidate if you've got had any of the subsequent conditions:

- Epilepsy
- Pregnant
- Breast Feeding
- Currently on light-sensitive medicine.
- Currently taking adrenal cortical steroid injections or the other quite steroid injection(s)
- LED light-weight medical aid is appropriate for those that would like to:
 - Minimize the looks of fine lines, wrinkles and crow's feet
 - Eliminate blemishes and minimize pores
 - Improve the looks of skin tone and texture
 - Restore the skin's natural albuminoidal production
 - Reduce animal pigment that causes brown spots
 - Reduce skin redness, scarring and dry patches caused by disease of the skin

1. Am i able to use disease of the skin merchandise with it?

Yes. Your results will certainly be quicker if you employ a top quality topical disease of the skin product. Associate disease of the skin product shouldn't be used before exploitation your light-emitting diode devices, however are often applied afterwards.

2. Am i intelligent to practice skin care stock in unification with it?

Yes. Light-weight medical assistance really wills intensification the uptake and use of your top eminence skin care commodities. Your results can happen quicker if charity using them.

3. However usually am i able to use it? And the way long can it desire see results?

- For the light-emitting diode devices, choose the color consistent with the condition you would like to treat. You'll be able to use it 20~25 minutes' sessions once on a daily basis, or 3~ four times per week. It depends on your age and current skin conditions.
- Most folks can see initial ends up in four ~8 weeks; several notice ends up in as short as one ~ three weeks; and a number of people can take as long as nine ~ eleven weeks.

I PLEASE continuously keep in mind THAT

1. Light-emitting diode light-weight technology works; it's simply a matter of "How long can it take". Everything takes time.
2. Even the foremost advanced system won't work if you aren't persistent, and this is often no exception. If you don't follow the directions and continue along with your treatments, you may not see the meant results and it should be straightforward to urge discouraged.

